

Transient Grating Spectrometer



Carrier diffusion coefficient
in a matter of minutes

Non-invasive
measurement technique

Fully automated and
computer controlled


Continuous setting of
grating period

Sensitivity down to $\mu\text{J}/\text{cm}^2$
excitation level

Advanced measurement
and analysis software

Photoluminescence (PL)
measurement option

HARPIA-TG principal scheme

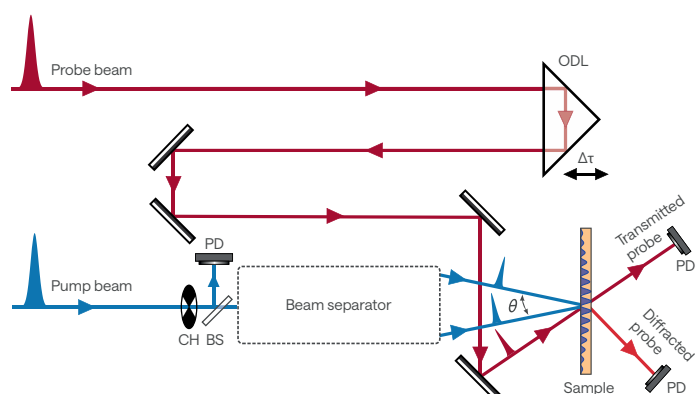

BS – beam splitter
CH – chopper
PD – photodiode

ODL – optical delay line
PM – parabolic mirror

θ – intersection angle
 $\Delta\tau$ – delay

Specifications

Measurement mode	Transmission	Reflection
Grating recording wavelength ¹⁾	340 – 560 nm	
Probe wavelength ²⁾	1030 nm	
Grating period ³⁾	1.05 – 12.5 μm	1.5 – 4.5 μm
Pulse duration	< 290 fs	
Delay range	Up to 8 ns	

MEASUREMENT RANGES

Diffusion coefficient	$\geq 0.1 \text{ cm}^2/\text{s}$
Carrier lifetime	3 ps – 8 ns


DIMENSIONS

Physical dimensions (L x W x H)	730 x 420 x 188 mm
---------------------------------	--------------------

¹⁾ Extendable to 750 nm by applying different physical gratings. Contact sales@lightcon.com for more details.


²⁾ OPA-based probe is available upon request. Contact sales@lightcon.com for more details.

³⁾ Depends on the pump wavelength.

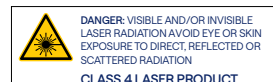


Performance

Diffusion coefficient of GaN
as a function of fluence

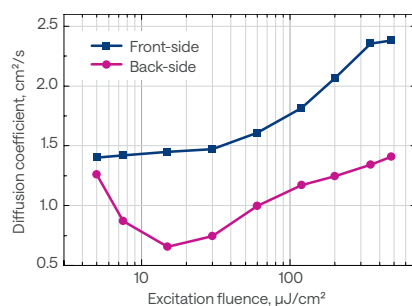

Carrier lifetime of GaN
as a function of fluence

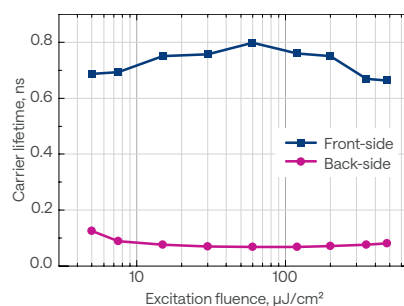

Diffusion length of GaN
as a function of fluence

