



MWM Laser Wavemeter



The Santec Laser Wavemeter is a precision compact laser wavelength measurement device, and optical spectrum analyser, with on-unit display and ethernet/USB connectivity standard. Use it standalone or with the included Windows GUI software package.

Unlike interferometric wavemeters, our device clearly reveals multimode laser operation, making it particularly suitable for use with external cavity diode lasers and atom cooling and trapping experiments. At prices so low you can afford to put a wavemeter on every laser in your lab, to know when your laser is unlocked or multimoding.

Features

- Picometre (<GHz) accuracy
- Picowatt sensitivity
- Any wavelength from 370 – 1560 nm[†]
- Instantly identifies multimode input
- Fast: up to 1250 measurements per second
- CW or pulsed laser input
- Built-in PID feedback with analogue output
- Ethernet and USB standard

Applications

- Laser frequency diagnosis and testing
- Identifying laser multimode behaviour
- Determining absorption lines
- Gas spectrometry
- Raman fluorescence

Laser Wavemeter

Specifications MWM

Wavelength/frequency

Wavelength range	370 – 1560 nm Supplied calibrated at one wavelength e.g. 780 ± 12 nm
Input power requirement	<1 pW (1 picowatt); max 30 nW
Precision	100 MHz (0.1 pm) at 780 nm
Accuracy	± 1 GHz / ± 0.001 nm at 780 nm, at time of calibration
Optical resolution	< 0.02 nm (σ std dev, wavelength dependent)
Dynamic range	> 35 dB (>50 dB with HDR)

Electronics

Display	Built-in LCD and host computer
Interface	Windows GUI
PID feedback	12-bit DAC output, 0.5 mV resolution, 10 Hz bandwidth
Readout	Typically 20 per second, up to 100 /s, host dependent

Interface

Ethernet	10/100 TP RJ45
USB	USB2.0, plug type USB-B (350 mA with display on)
SMA	Analogue output, ± 2.5 V, for PID feedback control
Power consumption	5 W


Inputs/outputs

Optical input	FC single mode fiber
DAC	12-bit output, ± 2.5 V, 0.5 mV resolution

Dimensions

Dimensions (approx.) 165 x 85 x 70 mm (LxWxH)

†Note: the MWM operates over a small wavelength range, e.g. 780 ± 12 nm. Operation in a different wavelength range requires mechanical adjustment and recalibration.

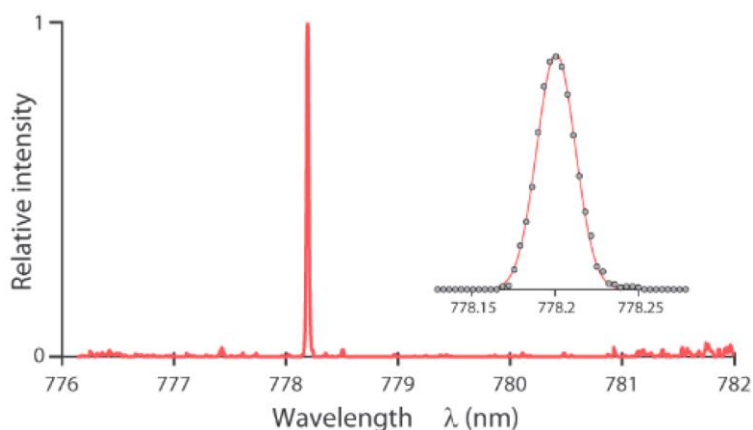

Photonic Solutions Ltd
Edinburgh
EH14 4AP
+44 131 664 8122
www.photonicsolutions.co.uk